



TURNER DESIGNS

845 W. Maude Avenue • Sunnyvale, CA 94085
(408) 749-0994 • FAX (408) 749-0998
www.turnerdesigns.com

*Picofluor*TM

Handheld Fluorometer



User's Manual

Dated: 03/19/2003

Version: 1.2

P/N 998-0853

Table of Contents

1. Introduction	
1.1 Description	4
1.2 Inspection and Setup	4
1.3 General Information and Precautions.....	5
2. Quick View Diagrams.....	6
3. Instrument Parameters	
3.1 Firmware Flowchart.....	7
3.2 Instrument Power Up	9
3.3 Fluorescence Channels.....	9
3.4 Calibration Standard Value	9
3.5 Calibration	10
3.6 Sample Analysis	10
3.7 Internal Data Logging (IDL)	11
3.7.1 Activate Data Logging	11
3.7.2 Download Data	11
3.7.3 Erase Data.....	11
3.8 Diagnostic Information	12
4. General Considerations for Analysis	
4.1 Handling Samples	13
4.2 Linear Range and Quenching	14
4.3 Temperature Considerations.....	15
4.4 Positioning Samples.....	15
4.5 Data Quality	15
5. Warranty	
5.1 Terms	16
5.2 Warranty Service	17
5.3 Out-of-Warranty Service.....	18

Appendices:

A. Instrument Specifications

A1. General Specifications.....	19
A2. Optical Specifications.....	19

B. Internal Data Logging

B1. Shipping Checklist	20
B2. Hardware Requirements	20
B3. Installation	20
B4. Connecting.....	20
B5. Troubleshooting	21

1. Introduction

1.1 Description

The *Picofluor*TM is a dual-channel fluorometer designed for quick, easy and accurate fluorescence measurements. When properly calibrated with a known standard, the *Picofluor*TM can be used to determine concentrations of unknown samples.

1.2 Inspection and Setup

1.2.1 Inspection

Upon receiving your instrument, please inspect everything carefully and make sure all accessories are present. All shipments include:

- The *Picofluor*TM
- The User's Manual
- 4 AAA batteries
- Internal Data Logging (IDL) package
 - PC interface cable
 - Spreadsheet Interface Software (1CD)

1.2.2 Setup

Before the *Picofluor*TM can be used, the supplied batteries must be installed.


1. On the back of the instrument, loosen the screw and remove the battery panel (see Section 2 for diagram).
2. Install the 4 AAA batteries into the appropriate spaces.
3. Replace the battery panel and tighten the screw. The panel has an o-ring, which creates a watertight seal. The battery panel may be difficult to install if there is no lubrication on

the o-ring. Use a silicon based grease to lightly lubricate the o-ring if necessary.


1.3 General Information and Precautions

- The sample compartment cannot accept glass or quartz cuvettes but readily accepts methacrylate or polystyrene cuvettes.
- A minimum volume of 2 mL in a 10 x 10 mm cuvette is required for best results.
- A minimum volume of 75 μ L in a Minicell adaptor is required.
- Avoid having any air bubbles in your sample. They can significantly affect the fluorescent reading.

2. Quick View Diagrams

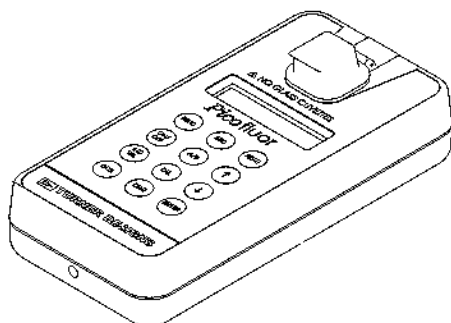



Battery Panel

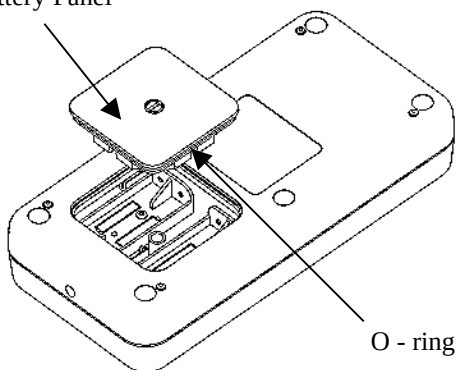

3. Instrument Parameters

3.1 Firmware Flowchart

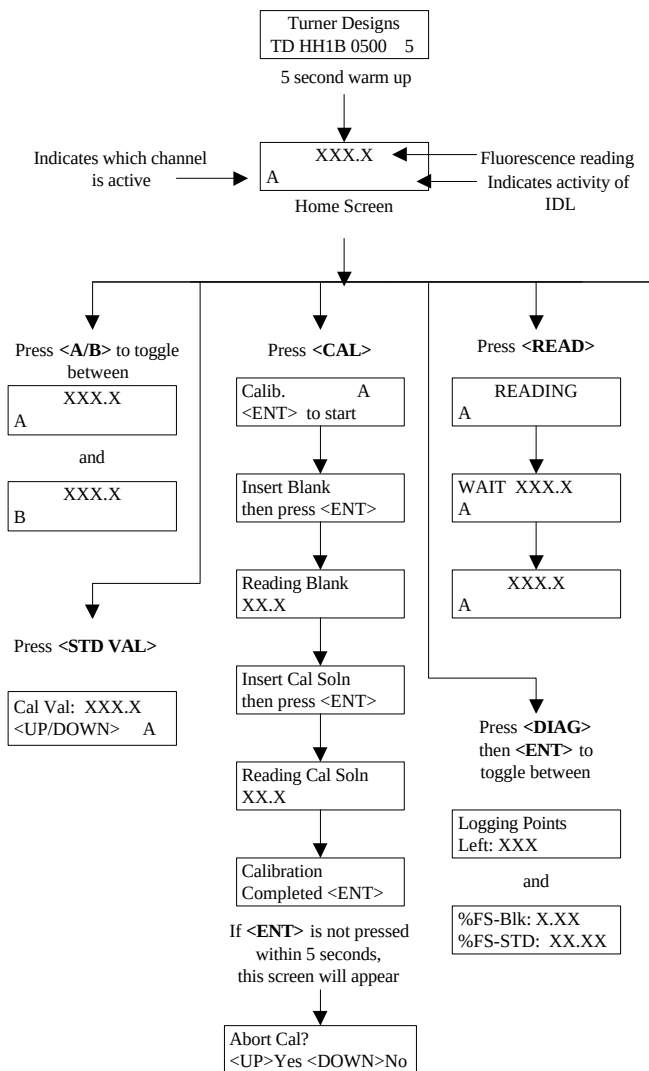

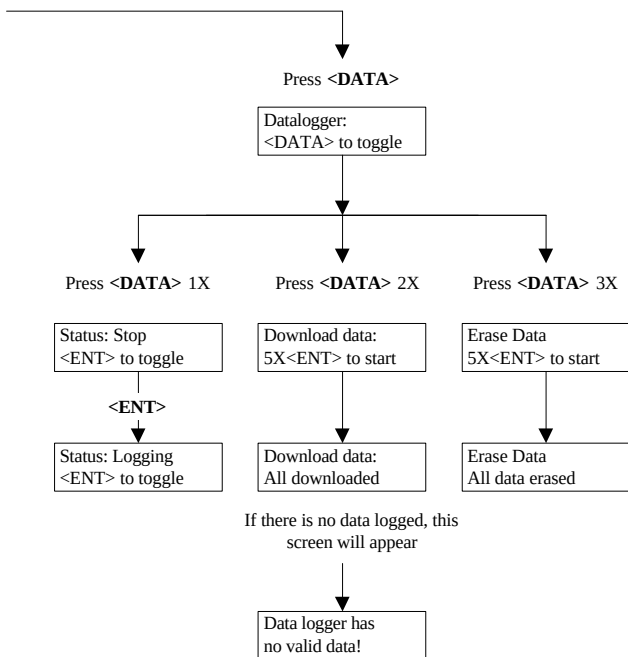

3.2 Instrument Power Up

To turn on the *Picofluor*[™], press the <ON/OFF> button. The instrument takes 5 seconds to warm up. After the warm-up, the *Picofluor*[™] is ready for operation.

3.3 Fluorescence Channels

Choose the appropriate channel for your analysis. To do this, press the <A/B> button to toggle between the 2 channels. The activated channel will be displayed in the lower left corner of the Home screen.

Model 8000-003

UV = Ultra Violet excitation
Blue = Blue LED excitation

Model 8000-004

RHOD = Green LED excitation
Blue = Blue LED excitation

See Appendix A for Optical (Wavelength) Specifications

3.4 Calibration Standard Value

Before performing a calibration, assign a numeric value between 1 and 999 for your standard solution.

1. Press the <STD VAL> button.
2. Use the up and down arrows to adjust the standard value. Holding this button down will activate faster scrolling.

3. When finished, Press <ESC> or <ENT> to accept the value and return to the Home screen.

3.5 Calibration

We recommend that you always calibrate before performing any sample analysis.

1. Press the <CAL> button.
2. Press <ENT> to start the calibration.
3. Insert your blank and press <ENT>. The *Picofluor*TM will average the blank fluorescence for 10 seconds.
4. Insert the calibration standard and press <ENT>. The *Picofluor*TM will average the standard fluorescence for 10 seconds.
5. Press <ENT> when the calibration is complete to accept the calibration. If <ENT> is not pressed within 10 seconds, you will be asked if you want to abort the calibration. Press the up or down arrow to abort or accept the calibration respectively.

If at anytime during steps 1-4 you want to stop the calibration, press <ESC>. This will return you to the Home screen and default back to the previous calibration.

3.6 Sample Analysis

1. Insert your sample.
2. Press either <READ> button. The instrument will autorange, then measure and average the fluorescence signal over a 5-second interval.
3. The result will be displayed at the top and center of the Home screen.

4. The top left corner will display “WAIT” for 5 seconds. Once “WAIT” disappears, another sample reading can be performed.

3.7 Internal Data Logging (IDL)

Your *Picofluor*[™] can log up to 1000 data points. The DATA screens control logging, downloading and erasing of data.

3.7.1 Activate Data Logging

1. Press the <DATA> button 2 times.
2. Press <ENT> to toggle between logging and stop.
3. Press <ESC> when finished to return to the Home screen.

3.7.2 Download Data

1. Connect the *Picofluor*[™] to the serial port of your PC using the provided interface cable.
2. Open the Spreadsheet Interface Software. (See Appendix B for computer requirements and installation).
3. Press the <DATA> button 3 times.
4. Press <ENT> 5 times to start the data download.
5. Press <ESC> when finished to return to the Home screen.

3.7.3 Erase Data

1. Press the <DATA> button 4 times.
2. Press <ENT> 5 times to erase all logged data.
3. Press <ESC> when finished to return to the Home screen.

3.8 Diagnostic Information

1. Press <DIAG> to access the diagnostic screens.
2. The first screen shows the number of data points available for internal data logging.
3. Press <ENT> to toggle to the %FS (Full Scale) values from the calibration blank and standard.
4. Press <ESC> when finished to return to the Home screen.


4. General Considerations for Analysis

4.1 Handling Samples

1. Take care not to spill samples into the sample chamber. Promptly wipe any spills.
2. The *Picofluor*TM is very sensitive and even small amounts of material from a previous sample may result in errors. Use a clean cuvette for each reading. Thorough and proper cleaning of cuvettes between sample readings is essential, and is especially important if you are using the same cuvette for samples and the blank.
3. Fill the cuvette at least 50% full; significant error will result if it is not filled enough.
4. The cuvette **MUST BE CLEAN AND DRY** on the outside when taking readings. Moisture and condensation on the outside can result in error.
5. Minute bubbles in samples will cause errors in the readings. Take care not to introduce bubbles into samples. Particular care must be taken with the Minicell adaptor. Slight tapping on the outside cuvette wall will often help dissipate bubbles.

4.2 Linear Range and Quenching

The linear range is the concentration range in which the readout of the *Picofluor*TM is directly proportional to the concentration of the fluorophore. The linear range begins with the smallest detectable concentration and spans to an upper limit (concentration) that is dependent upon the properties of the fluorescent material, the filters used, and the path length. See the diagram below.

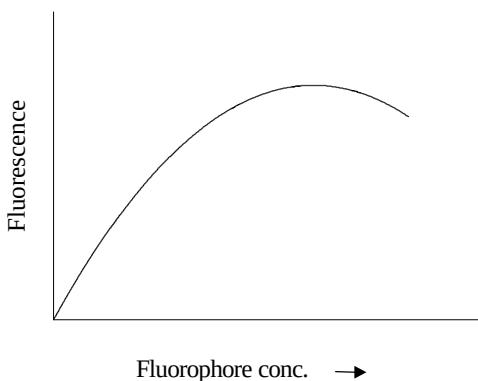

A nonlinear relationship is seen at higher concentrations where the fluorescence signal increases at decreasing rates in comparison to the change in concentration. At even higher concentrations, readings will begin to decrease even though the sample concentrations are increasing. This effect is known as "signal quenching."

Linearity may be checked by diluting a sample 1:1 or some other convenient ratio (be sure to use an appropriate blank for the dilutions). If the sample is within the linear range, the reading will decrease in direct proportion to the dilution. If the reading does not decrease in direct proportion to the dilution or if the reading increases, the sample is beyond the linear range of your fluorophore.

4.3 Temperature Considerations

Fluorescence is temperature sensitive. As the temperature of the sample increases, the fluorescence decreases. For greatest accuracy, read the blank, standard, and samples at the same temperature.

4.4 Positioning Samples

For low concentration samples, cuvettes often will give slightly different measurements depending upon their orientation in the sample compartment. This is due to defects in the shape of the cuvette that are not visible to the human eye. We recommend that the cuvette be marked at the top and positioned in the sample compartment the same way each time to minimize error.

4.5 Data Quality

The *Picofluor*[™] is only as accurate as the standards that are used to calibrate it. This is why it is important to take care when preparing standards, samples, and blank. One should follow good laboratory practices when preparing all solutions and reagents.

5. Warranty

5.1 Terms

Turner Designs warrants the *Picofluor*TM Fluorometer and accessories to be free from defects in materials and workmanship under normal use and service for a period of one year from the time of initial purchase, with the following restrictions:

1. The instrument and accessories must be installed, powered, and operated in compliance with the directions in this *Picofluor*TM User's Manual and directions accompanying the accessories.
2. Damage resulting from measurement of samples found to be incompatible with the materials used in the sample system is not covered.
3. Damage resulting from contact with corrosive materials or atmosphere is not covered.
4. Damage from salts, proteins and other moderately corrosive materials that are not promptly removed from the instrument are not covered.
5. Damage caused by modification of the instrument by the customer is not covered.

5.2 Warranty Service

To obtain service during the warranty period, the owner shall take the following steps:

1. Write or call the Turner Designs Customer Service or Applications Departments and describe as precisely as possible the nature of the problem.
2. Make minor adjustments or tests as suggested by Turner Designs personnel.
3. If proper performance is not obtained, ship the instrument, prepaid, to Turner Designs. The instrument will be repaired and returned free of charge for all customers in the contiguous continental United States.

For customers outside of the contiguous continental United States, and who have purchased our equipment from one of our authorized distributors, contact the distributor. If you have purchased direct, contact Turner Designs. We will repair the instrument at no charge, but we will not pay for shipment, documentation, etc. These charges will be billed at cost.

NOTE! Under no conditions should the instrument or accessories be returned without notice. Prior correspondence is needed:

- a. To ensure that the problem is not a trivial one, easily handled in your laboratory, with consequent savings to everyone.
- b. To specifically determine the nature of the problem, so that repair can be rapid, with particular attention paid to the defect you have noted.

5.3 Out-of-Warranty Service

Proceed exactly as for Warranty Service, above. If our service department can assist you by phone or correspondence, we will be glad to, at no charge.

Repair service will be billed on a basis of time and materials. A complete statement of time spent and materials used will be supplied. Shipment to Turner Designs should be prepaid. Your bill will include return shipment freight charges.

Address for Shipment:

Turner Designs
845 W. Maude Ave.
Sunnyvale, CA 94085

Appendix A: Instrument Specifications

A1. General Specifications

Specification	Description
Size	1.75" x 3.5" x 7.25" (4.45 cm x 8.9 cm x 18.4 cm)
Weight	14 oz (0.4 kg)
Dynamic range	4 orders of magnitude
Resolution	12 bits
LCD Display	2 x 16 characters
Case	Meets IP 67 Standard; dustproof and waterproof
Temperature	41—104 °F; 5—40 °C
Detectors	Photodiodes: measurement capability from 300—1000 nm
Calibration Type	Single-point and blank
Alarms	Low battery, circuit failure, high blank
Cuvette Type	10 x 10 mm plastic
Warm Up Time	5 seconds
Automatic Power Down	After 90 seconds of inactivity

A2. Optical Specifications

Light Source	UV LED	Blue LED	Green LED
Excitation Optics	365-395 nm	475±15 nm	525±20 nm
Emission Optics	>430 nm	515±20 nm	>570 nm
Limit of Detection		<50 pM Fluorescein	<200 pM Rhodamine

Appendix B: Internal Datalogging

B1. Shipping Checklist

Instruments with internal data logging will also receive in their shipment:

- Interface cable
- Spreadsheet Interface Software (1 CD)

Both of these items are necessary for retrieving the stored data in the *Picofluor*TM.

B2. Hardware Requirements

- PC with Windows 95 or later
- MS Excel 5.0 or later
- At least 1 available serial port

B3. Installation

1. Exit all Windows programs.
2. Insert CD and run the setup program.
3. The setup wizard will install the necessary files. When the setup is complete, an icon named "Spreadsheet Interface Software" will be found in the "Programs" menu.
4. Restart your computer.

B4. Connecting

1. Using the provided cable, connect the 9 pin adapter of the cable into the available serial port of your computer.
2. Plug the opposite end of the cable into the base of the *Picofluor*TM.
3. Open the Spreadsheet Interface Software.

4. Click on the box to the right of the COM port icon to select the appropriate COM port. This is usually COM port 1.
5. Click on “Start”. The program will find and open MS Excel and should be ready for data transfer. The boxes left of the COM port and MS Excel should both be green.
6. Follow the directions from Section 3.7 for collecting and downloading data from the *Picofluor*TM. Data will automatically appear in the excel spreadsheet. BE SURE to save this data BEFORE closing the software.

B5. Troubleshooting

Difficulties can arise when parameters are set incorrectly or connections with the cable are not tight enough. Here are some common problems.

1. Box to the left of the COM port is red. This means that the COM port is not available.
Causes:
 - a. Another instrument or program (such as palm pilot/ hot sync) could be occupying the port, making it unavailable. Make sure to close all programs of this type before downloading data
 - b. The port selected is incorrect. Follow step 4 of connecting to choose another COM port.
2. All lights are green, but no data transferred, even though the instrument says “All data downloaded”.
 - a. The connection between the instrument and the computer is bad. Check and tighten the cable connections. Make sure both ends of the cable are plugged in tightly.